


Abdominal compression bandage for horses

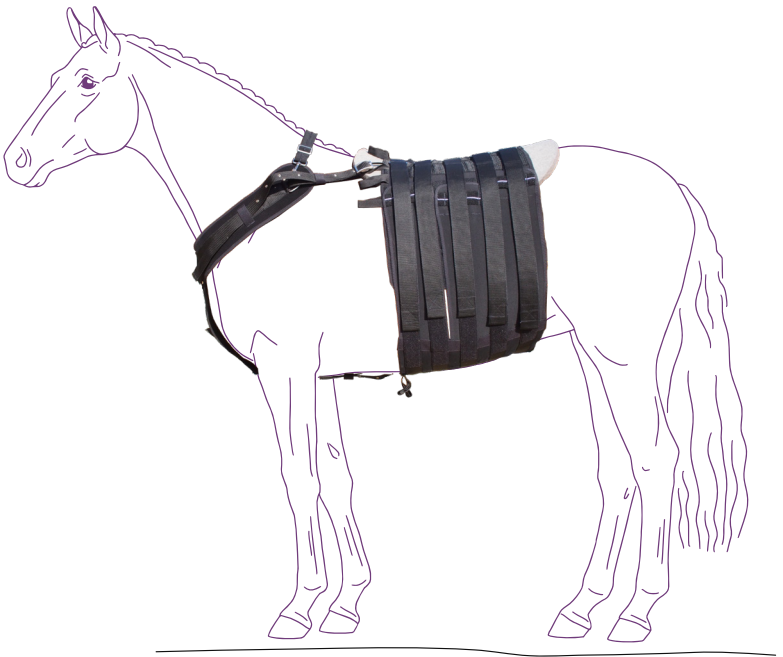

Can be fitted by one person
in just a few simple steps –
quick and very easy to adjust

The abdominal compression bandage from vetal

After colic surgery, the wound suture must be optimally protected and stabilized to prevent scar perforation, edema, and infection. However, a conventional bandage is particularly difficult to apply in this situation.

The **abdominal compression bandage** from vetal provides relief: It's quick and easy to apply, even by one person. Its fit is very flexible, fitting like a second skin, preventing common urination problems.

Abdominal compression bandage



This is achieved through a coordinated combination of partially elastic bandage material with firm stabilizing elements made of canvas (strong cloth).

The fastening is variable and can be done using Velcro straps, with which the compression can be continuously adjusted

The individually adjustable chest harness prevents the bandage from slipping backwards, its elastic parts maintain freedom of movement

To protect against sores, pressure and chafing, we also have a high-quality saddle pad developed. It is made of pure wool felt, which offers optimal moisture regulation and gently distributes pressure.



**stabilized
wound area**


The wound dressing material is placed in a replaceable cotton bag washable at 95°C. If additional compression is required, an additional compression plate can increase pressure in specific areas if indicated.



Facilitating daily control and supply: fixable wound dressing pocket



With medical recommendation:
Additional compression plate

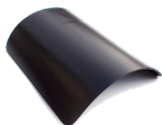

Our specially made saddle pad is designed for continuous use to protect the skin and coat.

For special back lines, they are available with a free spinal canal for the best comfort.



Abdominal compression bandage

Measurements & Order



.....
.....
Name & address of purchaser


.....
Phone & e-mail

Measurements

.....
Horse name/breed

.....
Height approx. age

.....
Diagnosis

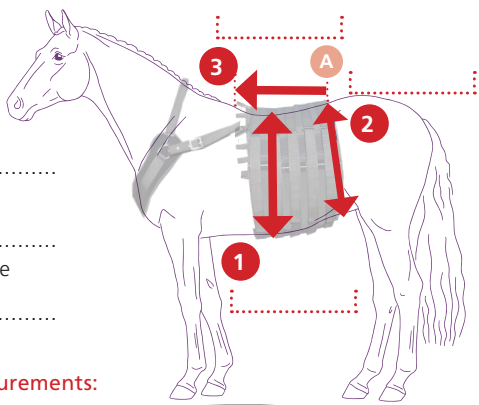

Please enter the 3 measurements:

- 1** measure the abdominal circumference at the widest point
- 2** measure the circumference at the narrowest point in front of the hindquarters. Point **A** where the measuring tape rests on the back,
- 3** is the starting point for the next measure.
From there, measure the length to the withers (= back measurement)

Order (please ask for the current price)

Amount

..... **Set**

Abdominal compression bandage +
chest harness + wound dressing pocket

.....
additionally

..... Saddle pad with free spinal canal

..... Compression plate

..... Wound dressing pockets (others for changing)

Send questions and your order to

vetal | Orthopedic technology for animals by albrecht GmbH
Chiemseestraße 81 | 83233 Bernau am Chiemsee | Germany
phone +49 8051 96129-0 | fax +49 8051 96129-30 | info@vet-al.de